

Distributed in New Zealand by:



Cooper Medical Limited
57 Mahunga Drive
Mangere Bridge
Auckland 2022



UROLOGY ISSUE 6
PRODUCT CATALOGUE

World Class Healthcare

Flexicare is a privately owned, leading manufacturer and supplier of medical devices with over 40 years of experience in the healthcare sector. Headquartered in the UK, Flexicare is a major supplier to the NHS and private institutions in the UK. We also have an international presence in over 105 countries through appointed distributors and supported by Flexicare subsidiaries in the US, Australia, China, Europe, India, Korea, Middle East, and Malaysia.



Intermittent catheters	4
Meatal dilators	7
Leg bags	8
F4 drainable 2 litre bag range	10
F4 Plus drainable bag	11
F2-EZ 2 litre bag range	12
F2 2 litre bag range	13
Urimeter	14
Wystopper stoma stopper	15
EZ-Flow catheter valve	16
FlexiSleeve leg bag holder	17
Bag stand	17
Leg bag straps	17
FlexiHang bag hanger	18

BLADDER DRAINAGE & STRICTURE THERAPY

Navigating even
the most difficult tract.



Flexicare has a wide range of intermittent catheters for bladder drainage and stricture therapy. We provide a range of Standard, Female and Paediatric lengths meeting all patients requirements.



Choice of Tips


Our selection of tips are available to navigate even the most difficult tracts and are available in Standard, Tiemann, Tapered and Soft Tip.

Hydrophilic or Uncoated

Our sterile catheters are supplied in a choice of hydrophilic or uncoated catheters in easy to open packaging.

Colour Coded

All our catheters and meatal dilators are colour coded for easy size recognition.

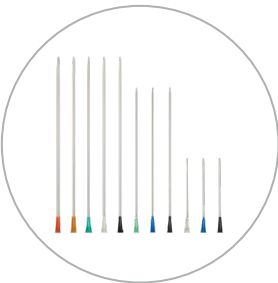

Step-by-Step Catheter Guides

Our easy to follow adult and baby guides help patients keep track of catheterisation and guide them through the procedure with ease. Visit our download and video pages at myflexicare.com.

WYCATH HYDROPHILIC INTERMITTENT CATHETERS

Our sterile Hydrophilic coated catheter comes with our Wycath water sachet sealed inside the packaging. Simply press the button to release the water before insertion. Easy to open for patients with dexterity difficulties and supplied in a range of tips.

Standard Hydrophilic Intermittent Catheters

TIP	LENGTH	EYELETS	● 10CH	○ 12CH	● 14CH	● 16CH	● 18CH
Soft Tip	40cm	2	W2SM10	W2SM12	W2SM14	W2SM16	W2SM18
Soft Tip	40cm	4	W4SM10	W4SM12	W4SM14	W4SM16	-
Tapered	40cm	2	W2TAM10	W2TAM12	W2TAM14	W2TAM16	W2TAM18
Standard	40cm	2	W2M10	W2M12	W2M14	W2M16	-
Tiemann	40cm	2	W2TM10	W2TM12	W2TM14	W2TM16	-

Female Hydrophilic Intermittent Catheters

TIP	LENGTH	EYELETS	● 8CH	● 10CH	○ 12CH	● 14CH
Soft Tip	18cm	2	W2SF08	W2SF10	W2SF12	W2SF14
Tapered	18cm	2	W2TAF08	W2TAF10	W2TAF12	W2TAF14
Standard	18cm	2	W2F08	W2F10	W2F12	W2F14
Standard	11cm	2	W2SH08	W2SH10	W2SH12	W2SH14

Paediatric Hydrophilic Intermittent Catheters

TIP	LENGTH	EYELETS	● 6CH	● 8CH	● 10CH
Soft Tip	30cm	2	-	W2SP08	W2SP10
Tapered	30cm	2	W2TAP06	W2TAP08	W2TAP10
Standard	30cm	2	W2P06	W2P08	W2P10
Tiemann	30cm	2	W2TP06	W2TP08	W2TP10



WYCATH UNCOATED INTERMITTENT CATHETERS

Our extensive choice of tips are available across the WyCath range, providing the ideal catheter for each patient.

Wycath’s sterile uncoated intermittent catheters are supplied in Standard, Female and Paediatric lengths meeting all patient’s requirements.

Standard Uncoated Intermittent Catheters

TIP	LENGTH	EYELETS	● 10CH	○ 12CH	● 14CH	● 16CH
Soft Tip	40cm	2	USM10	USM12	USM14	USM16
Soft Tip	40cm	4	U4SM10	U4SM12	U4SM14	U4SM16
Tapered	40cm	2	UTAM10	UTAM12	UTAM14	UTAM16

Female Uncoated Intermittent Catheters

TIP	LENGTH	EYELETS	● 8CH	● 10CH	○ 12CH
Soft Tip	18cm	2	USF08	USF10	USF12
Tapered	18cm	2	UTAF08	UTAF10	UTAF12

Paediatric Uncoated Intermittent Catheters

TIP	LENGTH	EYELETS	● 5CH	● 6CH	● 8CH	● 10CH
Soft Tip	30cm	2	-	-	USP08	USP10
Tapered	30cm	2	-	UTAP06	UTAP08	UTAP10
Standard	30cm	2	WUP05*	UP06	UP08	-
Tiemann	30cm	2	-	UTP06	UTP08	-

*This catheter is 25cm in length.



Standard Tip



Tiemann Tip



Tapered Tip



Soft Tip



Our market-leading meatal dilators are designed to ease through strictures. Available in a choice of 2 tips, with a hydrophilic coating to ease insertion.



Choice of Tips

Our meatal dilators are available in Standard and Tapered tip to help ease through strictures.

Hydrophilic Coating

Our Meatal Dilators have a hydrophilic coating to ease insertion. Simply fill the dilator sleeve with water then wait 30 seconds for hydrophilic coating to activate.

Meatal Dilators - 11cm Length

TIP	● 8CH	● 10CH	○ 12CH	● 14CH	● 16CH	● 18CH	● 20CH
Tapered	WMT08	WMT10	WMT12	WMT14	WMT16	WMT18	WMT20
Standard	WCM08	WCM10	WCM12	WCM14	WCM16	WCM18	-

CONTOURED, DISCREET LEG BAGS

offering a secure
and comfortable fit.



Our range of leg bags are contoured to offer a comfortable fit whilst distributing urine evenly, giving patients confidence wherever they are.

Non-return Valve

All our drainage bags are fitted with a non-return valve (NRV) as standard.

Date and Use Box

Clear indications of use and a date box allow information to be recorded on all our bags.

Range of Lengths

Available in 350ml, 500ml and 750ml with a choice of tube lengths.

Needle-Free Sample Port

Our sample port eliminates the risk of needle stick injury and the smooth exterior reduces potential marking of the skin and visibility under clothing.

Secure Connection

Our stepped connectors are featured on all our leg bags offering secure and reliable connection to catheters and drainage bags.

Easy to Operate, Taps

The lever taps fitted to our leg bags are 180 degree, 2 stage action taps. Easy to operate and reducing the risk of accidental opening.



LEG BAGS

Our leg bags are available in short (7cm) or long (30cm) tube lengths and for those requiring a customised length, we offer a 500ml capacity bag with an extra-long adjustable length tube (60cm).

Standard Leg Bags

Our standard leg bag is contoured offering an improved and comfortable fit while distributing urine evenly.

DESCRIPTION	SHORT TUBE HOSPITAL	LONG TUBE HOSPITAL	SHORT TUBE COMMUNITY	LONG TUBE COMMUNITY	ADJUSTABLE COMMUNITY
350ml Sterile leg bag	-	-	00-1352	00-2352	-
500ml Sterile leg bag	00-1503L	00-2503L	00-1502	00-2502	00-3502
500ml Non-Sterile leg bag	-	-	-	-	00-3501
500ml Sterile leg bag + gloves	-	-	00-1502G	00-2502G	
750ml Sterile leg bag	-	-	00-1752	00-2752	-

Discreet Leg Bags

Our Discreet leg bags are a neutral colour which makes them less conspicuous under clothing and the white print reduces the possibility of bright colours showing through clothing.

DESCRIPTION	SHORT TUBE COMMUNITY	LONG TUBE COMMUNITY	ADJUSTABLE COMMUNITY
350ml Sterile leg bag	00-5352	00-6352	-
500ml Sterile leg bag	00-5502	00-6502	00-6504
500ml Non-Sterile leg bag	-	-	00-6505
500ml Sterile leg bag + gloves	00-5502G	00-6502G	-
750ml Sterile leg bag	00-5752	00-6752	-

Hospital codes with non-flocked backing are supplied in packs of 100 and community codes with flocked backing for improved comfort in packs of 10.

giving patients
complete confidence
wherever they are



DRAINABLE OVERNIGHT BAGS

robust and versatile
drainage solution.



Our drainable urine bags provide a robust and versatile urine drainage solution for both hospital and community environments. All drainable 2 litre bags are fitted with a non-return valve to prevent urine reflux.

Easy to Operate, Taps

The F4 bags are fitted with a choice of a lever tap or T-tap. Both taps are simple and easy to operate.



Free and Rapid Flow

Our soft wide bore drainage tubing allows unimpaired flow of urine and comfort against the skin.

Lower Profile

The F4LS side tap lowers the profile of the drainage bag and moves the tap clear off the floor to prevent accidental contamination.

Needle-Free Sample Port

The needle-free sample port eliminates the risk of needle stick injury.

DESCRIPTION	HOSPITAL	COMMUNITY
F4 sterile 2 litre drainable bag with T-tap	00-1204	00-1200
F4L sterile 2 litre drainable bag with lever tap	00-1204L	00-1201
F4LS sterile 2 litre drainable bag with lever tap and side pouch	00-1207L	-
F3 non-sterile 2 litre drainable bag with T-Tap	00-2206	-

Hospital codes are supplied in packs of 100 and community in packs of 10.

F4 PLUS 2 LITRE BAG RANGE

The F4 Plus bag is primarily used with patients who have undergone surgery. It has an integral hanger and a pasteur drip chamber to observe urine drainage.

Hydrophobic Filter

The F4 Plus hydrophobic filter prevents air building up in the bag without allowing proteins or fluids to pass through.

Integral Bag Hanger

The hanger can fit securely to any bed. The inlet tube is maintained in position to allow for direct flow into the drip chamber.

Non-Return Valve

The F4 Plus bag is fitted with a second non-return valve preventing the contents of the bag returning to the drip chamber, eliminating reflux and providing directional flow of fluids.

Drip Chamber

The pasteur drip chamber enables the urine flow to be observed.

Enclosed Tap Protector

The lever tap can be placed in the protective tap holder which keeps the outlet clear of possible contamination from the floor whilst protecting against accidental opening.



DESCRIPTION	HOSPITAL	QTY
F4 Plus sterile 2 litre drainage bag	00-1208	20



F2-EZ is a 2L drainage bag designed to be opened quickly whilst allowing controlled drainage of fluids. The bag is designed with a simple twist-snap tap to enable the urine to be directed and emptied easily.



F2-EZ Top Located Tap

The Twist-snap tap is located in the top corner of the bag to offer a 'jug' effect and greater control when emptying.



Secure Connection

Our stepped connectors are featured on the F2-EZ and F2-EZB offering secure and reliable connection to catheters and drainage bags.



Twist-Snap Tap

The innovative easy-to-open angled connector with wide bore for rapid drainage is twisted to break the seal and drain.

DESCRIPTION	HOSPITAL	COMMUNITY
F2-EZ non-sterile 2 litre once drainable bag with top tap	00-2203	00-2203C
F2-EZB non-sterile 2 litre once drainable bag with bottom tap	00-2208	00-2208C

Hospital codes are in packs of 100 and community in packs of 10.

F2 2 LITRE BAG RANGE

Flexicare F2 bag provides a clean and safe method of urine drainage. F2 is a high quality bag that is simple and easy to use, assuring confidence to clinicians and comfort for users.

Easy to use

The simple and easy-to-use tear mark can be used to open and drain the bag when needed. After use, the bag can be easily disposed of making the F2 an ideal bag for use in both hospital and community settings.

Safe and Secure Drainage

The non-return valve (NRV) is designed to allow low pressure drainage in one direction while preventing all flow in the opposite direction, reducing the risk of ascending infection.

Secure Connection

The stepped connector is easily inserted into your drainage system and provides a secure connection.



DESCRIPTION	HOSPITAL	COMMUNITY
F2 non-sterile 2 litre once drainable tear to open bag	00-2201*	00-2202C
F2 sterile 2 litre once drainable tear to open bag	00-2202	-

*Non-sterile hospital codes are supplied in packs of 250, sterile codes are in packs of 100 and community in packs of 10.

easy to use
simply tear the bag
drain and dispose

URIMETER

The Flexicare Urimeter is a sterile closed collection system with an integral measuring chamber. The vertical drainage system and incorporated non return valves (NRV's) ensure the measuring chamber empties fully once the reading has been taken to ensure accurate measurements especially with low urine output.



Enclosed Tap Protector

The lever tap can be housed in the holder when not in use to prevent the tap touching the floor or being accidentally opened.

Hydrophobic Filters

The hydrophobic filter allows for the air and any gas build up in the bag to escape without allowing proteins or fluids to pass through.

Safe and Secure Drainage

Non-return valves reduce the risk of ascending infection and allows low pressure drainage without compromising the anti-flux property.

Clear print and date box

Clear indications of use and a date box allow information to be recorded and use monitored.

Vertical Drainage

The vertical drainage ensures accurate measurements.

Secure fastening

The integral hangers and straps ensure the Urimeter is securely attached to any bed frame.

DESCRIPTION	HOSPITAL	QTY
Urimeter with 150mm Tubing	00-1227	20



Wystopper- a flexible clear soft stoma stopper designed by a clinician to maintain stoma sites, especially continent urinary diversions and antegrade continence enema (ACE) tracts.




Tapered Tip

The clear soft silicone aids patient comfort and tapered tip to ease insertion.



Mushroom Cap

The mushroom shaped cap lies flush against the skin for discreetness.

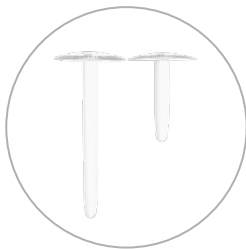

Choice of sizes

Available in sizes 10Ch, 12Ch, 14Ch, and 16Ch in 30mm or 60mm lengths.

DESCRIPTION	10Ch	12Ch	14Ch	16Ch	QTY
Wystopper Stoma Stopper 30mm	WS0310	WS0312	WS0314	WS0316	3
Wystopper Stoma Stopper 60mm	WS0610	WS0612	WS0614	WS0616	3

EZ-FLOW CATHETER VALVE

discreet and secure



EZ-Flow is an inconspicuous light weight neutral coloured catheter valve. The valve is designed to offer patients the flexibility of choice in managing their urinary drainage system. In line with all catheter valves, the valve remains attached to the catheter for 5 – 7 days.

Patient Comfort

The outlet sleeve is made of soft PVC that reduces pressure marks or discomfort for the patient.

Secure Fitting

The Ez-Flow multiple ridged connector ensures a secure fit with the catheter.

Patient Suitability

Suitable for patients with adequate bladder capacity, sensation, good dexterity and cognitive awareness.

Easy to Operate, Lever Tap

Ez-Flow has a 180°, 2 stage action tap. Easy to operate, the slightly firmer second stage reduces the risk of accidental opening on clothing or bedding.



DESCRIPTION	CODE	QTY
EZ-Flow Catheter Valve	00-0060	5

FLEXISLEEVE LEG BAG HOLDER

FlexiSleeve is a knitted Lycra, latex free leg bag holder to replace leg straps.

Distributing weight evenly

Offering more comfort to the bag wearer by distributing weight evenly around the leg.



DESCRIPTION	SMALL	MEDIUM	LARGE	X-LARGE	QTY
Leg Bag Holder	00-0035	00-0036	00-0037	00-0038	3


BAG STAND

Designed to hold any 2 litre bag ensuring the tap is clear off the floor to reduce the risk of potential contamination.

Plastic Coating

The strong metal stand is plastic coated making it easy to clean and ensuring that there are no sharp edges.

DESCRIPTION	CODE	QTY
Bedside bag stand	00-0022	50

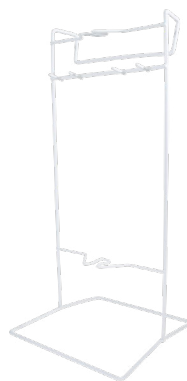

LEG BAG STRAPS

Flexicare leg straps feature an anti-slip silicone weave which grips the leg and ensures a secure and comfortable fit.

Secure Fastening

Secure Velcro fastening that can be easily adjusted.

DESCRIPTION	HOSPITAL	COMMUNITY
Lycra Anti-Slip Leg Straps	00-0032	00-0032C
Neutral Lycra Anti-Slip Leg Straps	-	00-0033C



Hospital codes are supplied in packs of 100 and community in packs of 10.

FLEXIHANG BAG HANGER

secure overnight
drainage



FlexiHang bag hanger is a universal 2 litre urine bag hanger uniquely designed to be fitted to beds with or without side bars. FlexiHang can be secured between the mattress and bed frame or hung from the bed using the flexible hooks.

Flexible Hooks

Wide flexible hooks ensure stability whilst hooked onto the bed.

Tube Holder

Secures the outlet tube in a vertical position to ensure optimum urine flow.

Easy Clips

The two snap-fit clips shut tight to hold the urine bag firmly in place.




DESCRIPTION	HOSPITAL	COMMUNITY
FlexiHang Bag Hanger	00-0014	00-0012

Hospital codes are supplied in packs of 100 and community in packs of 10.



Our website
myflexicare.com

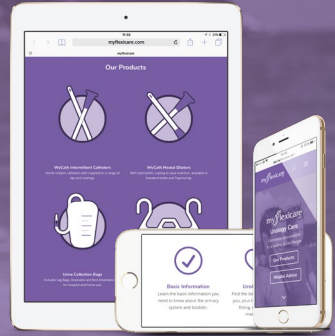

Built with you in mind.

Our website has been designed to offer the most recent information on our product range and new developments, providing you with all the information you need in one location.

Product guides and videos

Our easy to follow step by step guides are full of information on your leg bag, catheter valve and drainage bag. Download our personal guides or watch our videos designed to show you how to use our products with ease.

Helpful advice and tips

Discover our information and helpful advice on your catheter, leg bag or other continence care product and how to live a healthy, active lifestyle alongside your urology product.



Download your
personal guides



Watch videos of
products in use



Helpful tips for an
active lifestyle



Request your
samples online



Flexicare (Group) Limited

Cynon Valley Business Park, Mountain Ash, CF45 4ER

☎ +44 (0)1443 474647 ☎ +44 (0)1443 474222 ✉ enquiries@flexicare.com

Flexicare, Inc.

Irvine | CA | USA

☎ +1 949-450 9999

✉ inquiries@flexicare.com

Flexicare Middle East & Africa

Amman | Jordan

☎ +962 6554 8944

✉ enquiries@flexicare.com

Flexicare Australia Pty Ltd

Victoria | Australia

☎ +61 450 977 007

✉ enquiries@flexicare.com

Flexicare Medisize Germany

Siegburg | Germany

☎ +49 2241 93 86 0

✉ info@medisize.de

Flexicare Medisize Netherlands

Hoofddorp | The Netherlands

☎ +31 (0)88 77 44 555

✉ info@medisize.com

Flexicare Medisize Italy

Rome | Italy

☎ +39 065935621

✉ it_info@medisize.com

Flexicare Medical India Pvt. Ltd

Chennai | India

☎ +91 44 4554 0220

✉ enquiries@flexicare.com

Flexicare Korea Co., Ltd

Seoul | Republic of Korea

☎ +82 2 6949 3622

✉ info.kr@flexicare.com

Flexicare Medical Dongguan Limited

Dongguan City | China

☎ +86 769 83713371

✉ enquiries@flexicareasia.com

Prisma Flexicare Sdn Bhd

Kuala Lumpur | Malaysia

☎ +603-2280 6948

✉ enquiries@flexicare.com

Flexicare Lanka Pvt Ltd

Colombo | Sri Lanka

☎ +61450977007

✉ enquiries@flexicare.com

PT Flexicare Medical Indonesia

Jakarta | Indonesia

☎ +6010 222 9161

✉ enquiries@flexicare.com

www.myflexicare.com

FC4846 | VP5500-58B | 102019

Ph: 0800 800 823 - www.coopermedical.co.nz - admin@coopermedical.co.nz